



Android Settings for EasySense, EasyTouch and EasyPatch (Redmi 10X)

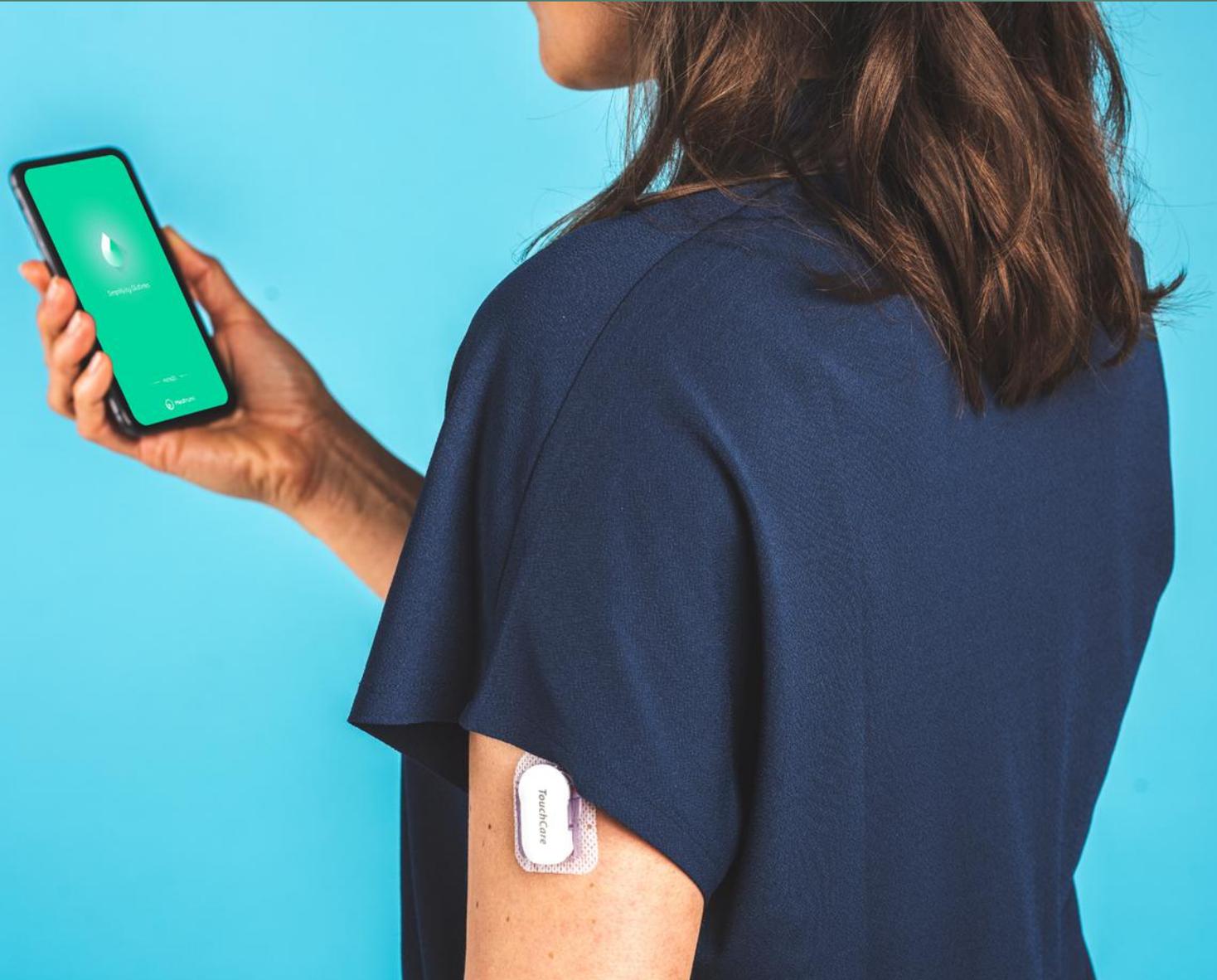

1. Connections Settings

- Make sure your smart device's **Bluetooth** is available and turn it on.

Settings→**Bluetooth**

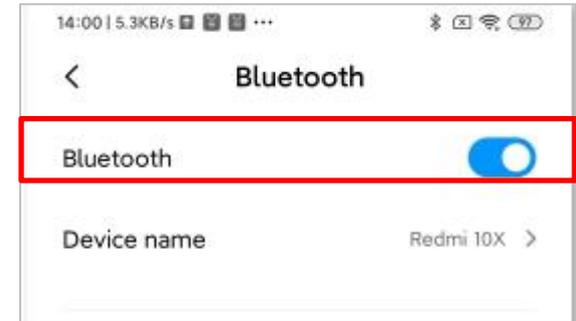

- Connect your smart device to Internet via WLAN or mobile data.

Settings→**WLAN**





2. Notifications Setting

- Turn on **Notifications** and allow Medtrum App to show notifications on the status bar.
- Turn off **Silent** and **Do Not Disturb** functions

Settings → Notifications → Medtrum apps

or

Settings → Apps → Medtrum App → Notifications

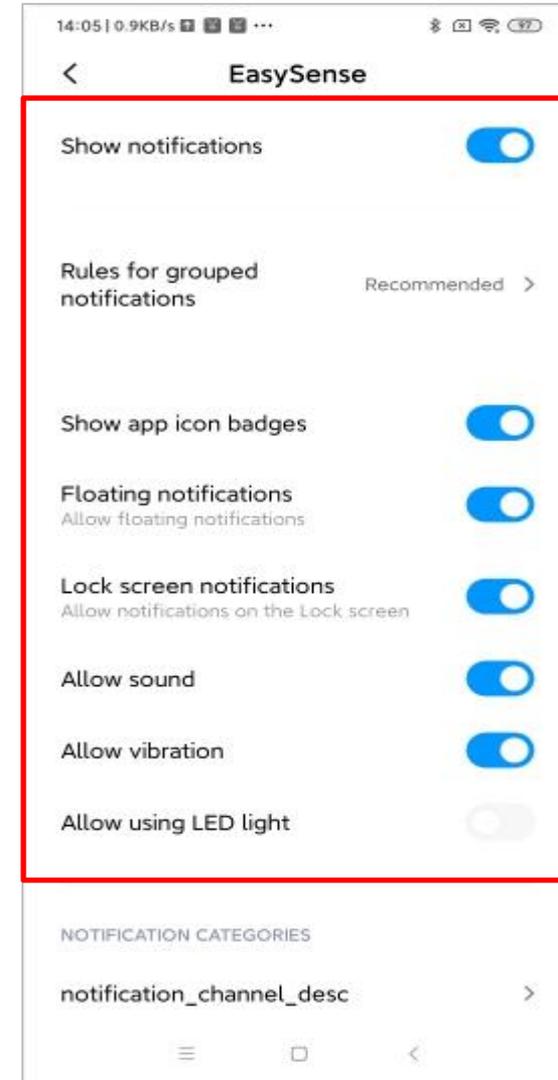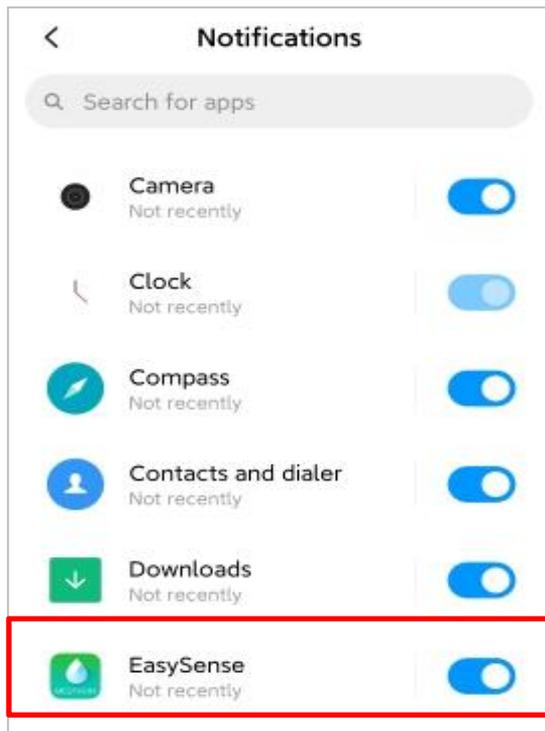


3. Location Setting

For Android system 10.0 and 11.0:

- Turn on **Location** so that the Bluetooth feature can be used.

Settings→**passwords&security**→**Location**



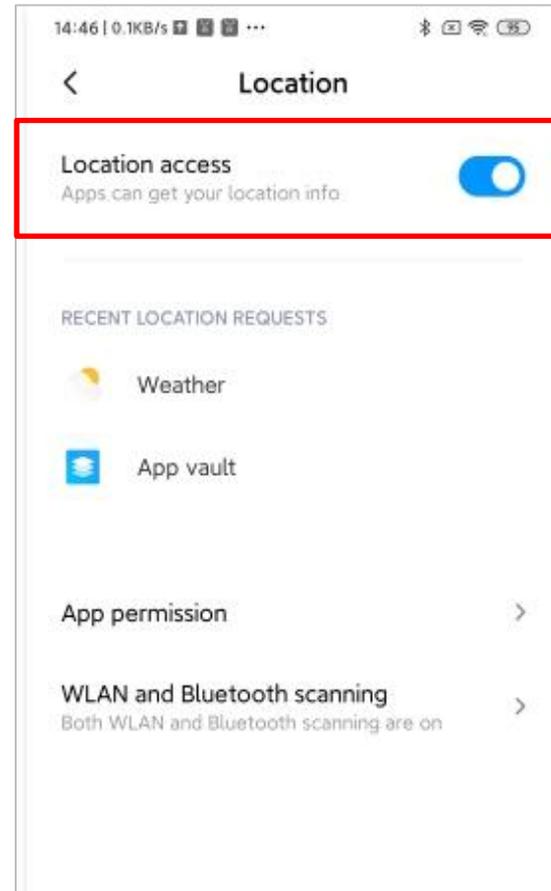


4. App Permissions Setting

For Android system 10.0 and 11.0:

- Allow Medtrum App to use **Camera, Location and Storage**

Settings→**Apps**→**ManageAPP**→**Medtrum App**→**Permission**

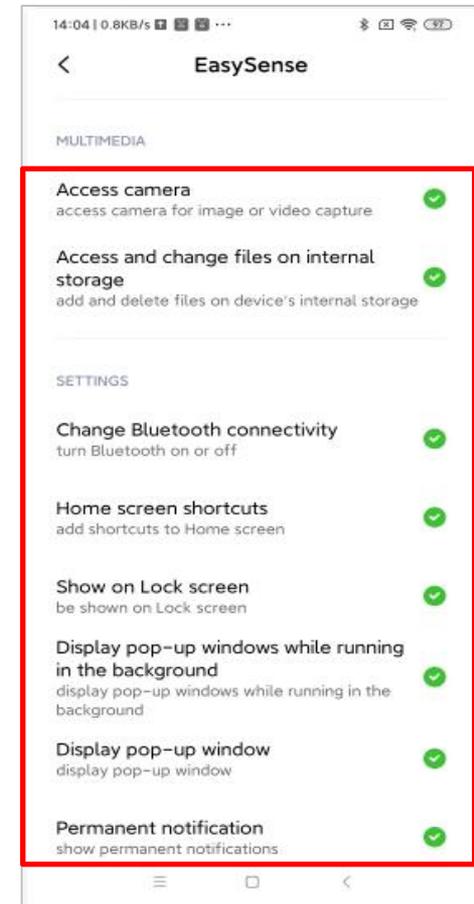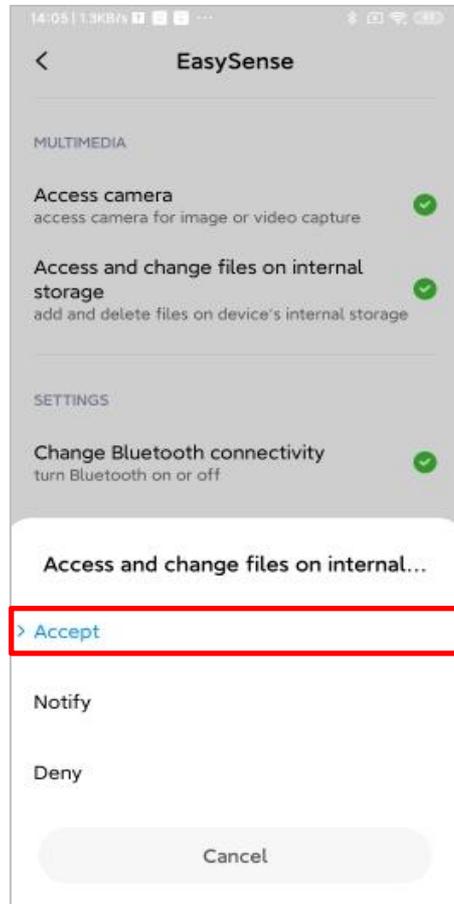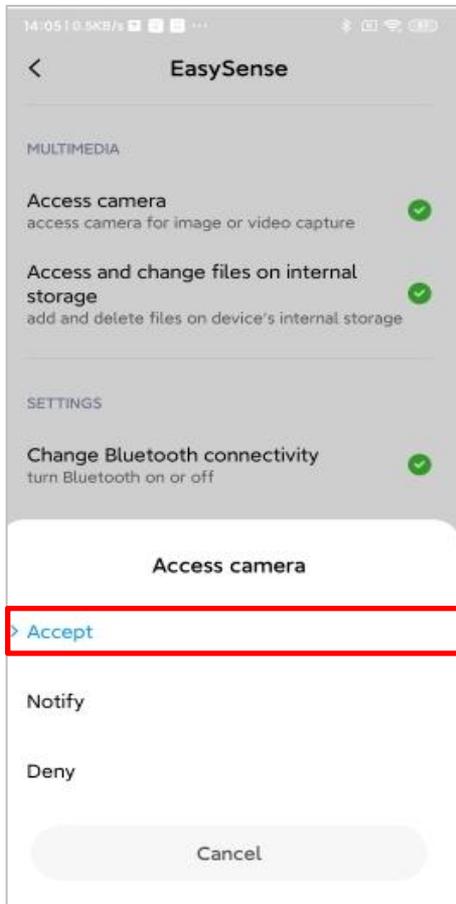


4. App Permissions Setting

For Android system 10.0 and 11.0:

Settings→**Apps**→**Medtrum App**→**Permissions**




Note: Allow **Location** to be used all the time.

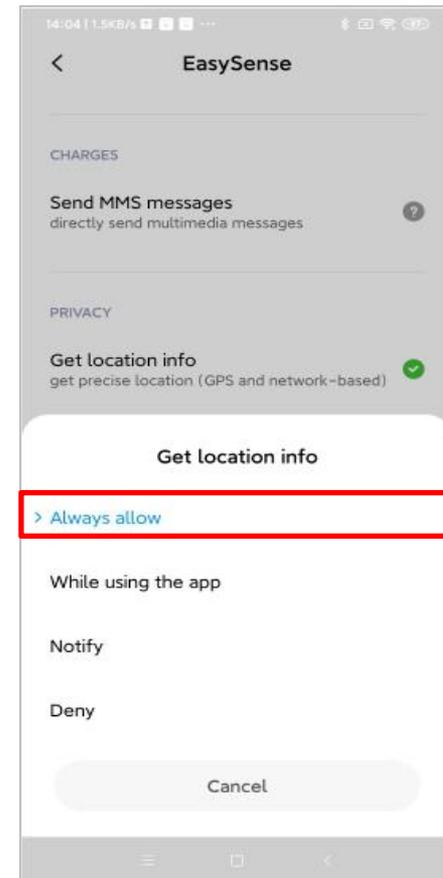

4. App Permissions Setting

For Android system 12.0 and above:

- Allow Medtrum App to use **Camera, Nearby Devices** and **Storage**
Settings→**Apps**→**ManageAPP**→**Medtrum App**→**Permission**

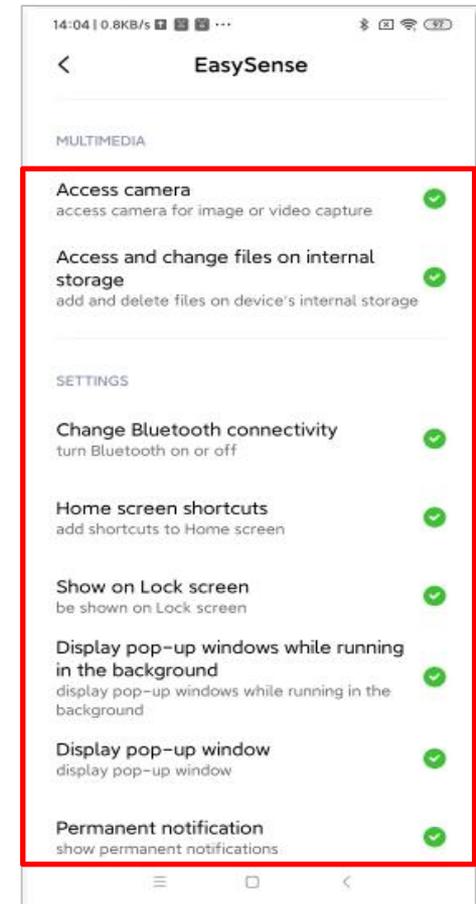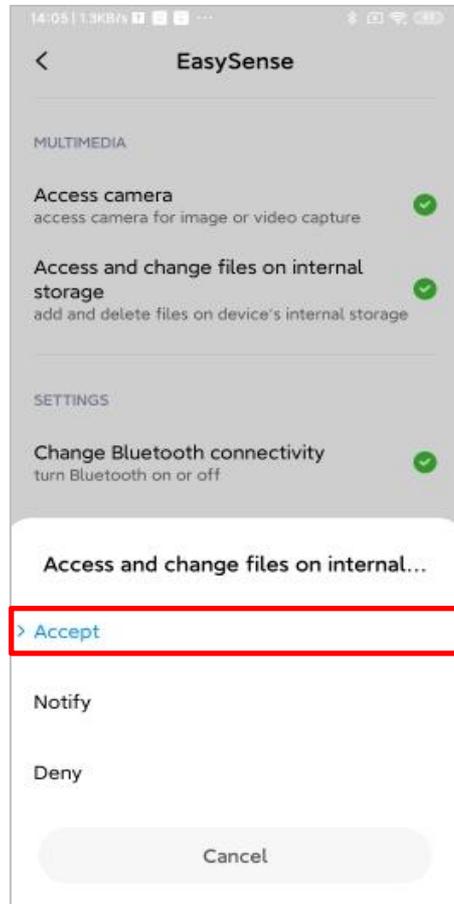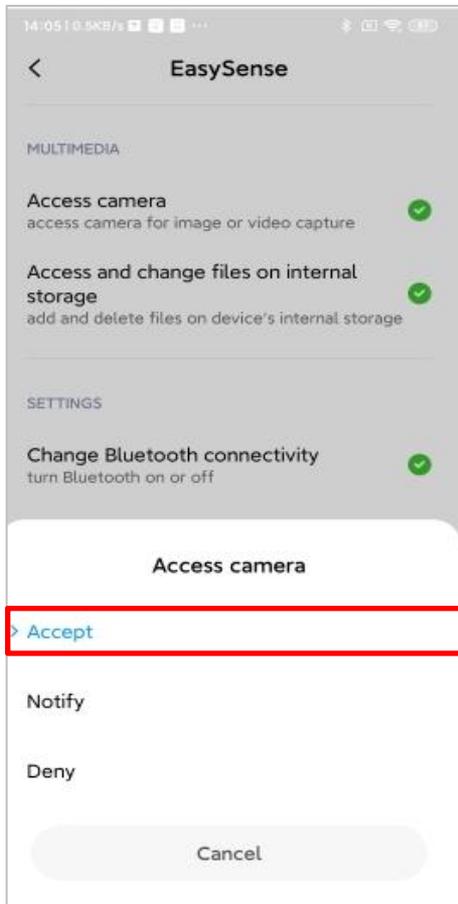


4. App Permissions Setting

For Android system 12.0 and above:


Settings→**Apps**→**Medtrum App**→**Permissions**

Note: Allow **Nearby Devices** to be used all the time.

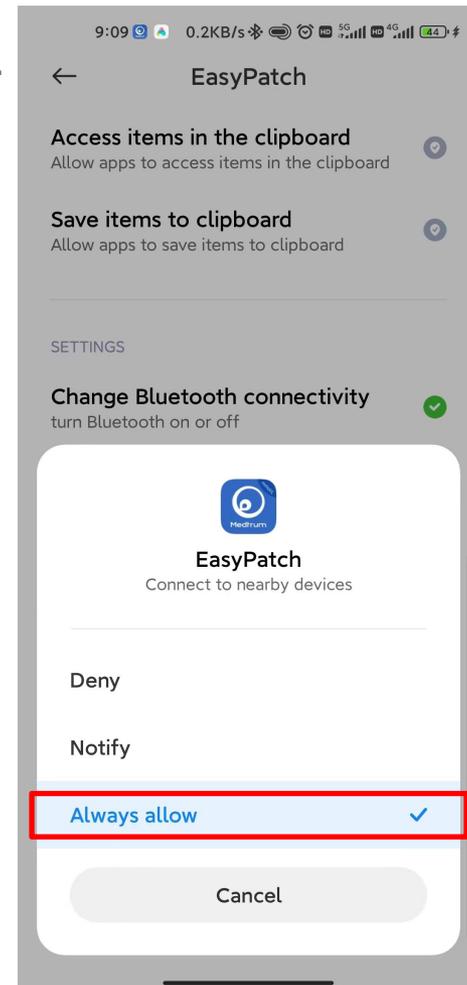

5. Battery Settings



- Turn off **Battery Saver** and **Ultra battery Saver**.
Settings→**Battery & Performance**

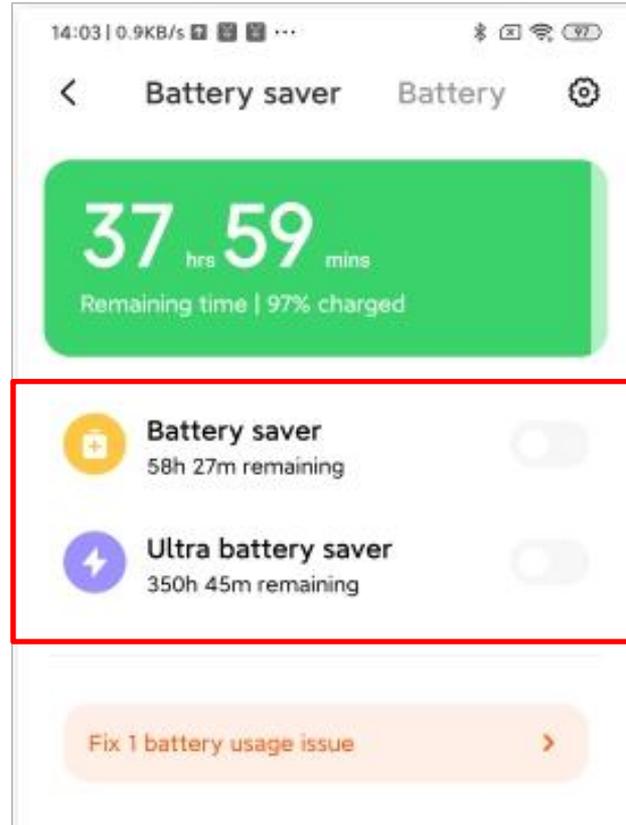

5. Battery Settings

- Select **Never** for **Turn mobile data when device is locked** .
- Select **Never** for **Clear cache when device is locked**.


Settings→**Battery & Performance**→**Setting icon**

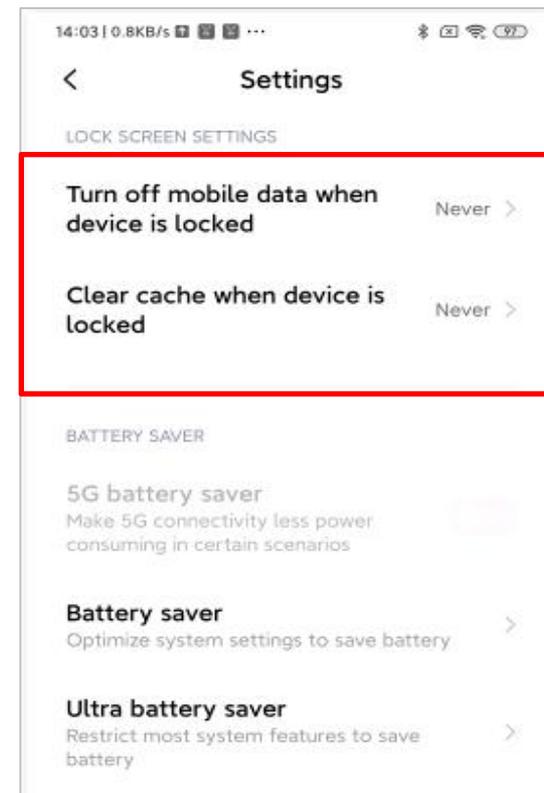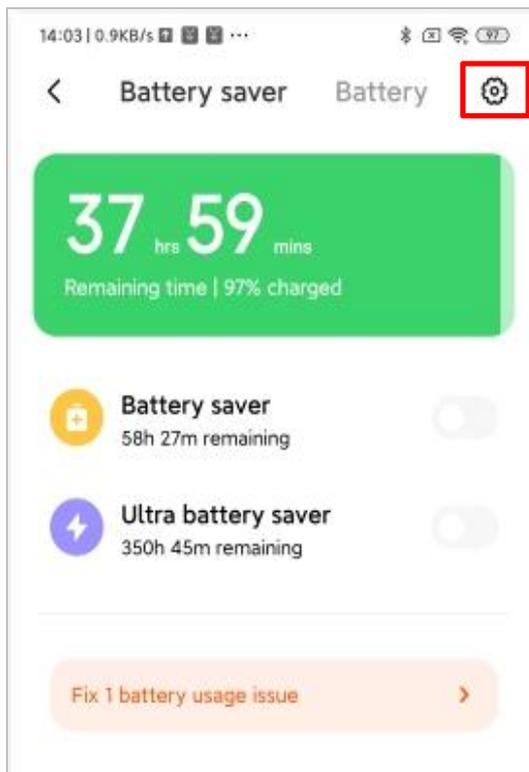

5. Battery Settings

- Select **No restrictions** for **Background settings**


Settings→**Apps**→**Manage apps** →**Medtrum App**→**Battery saver**

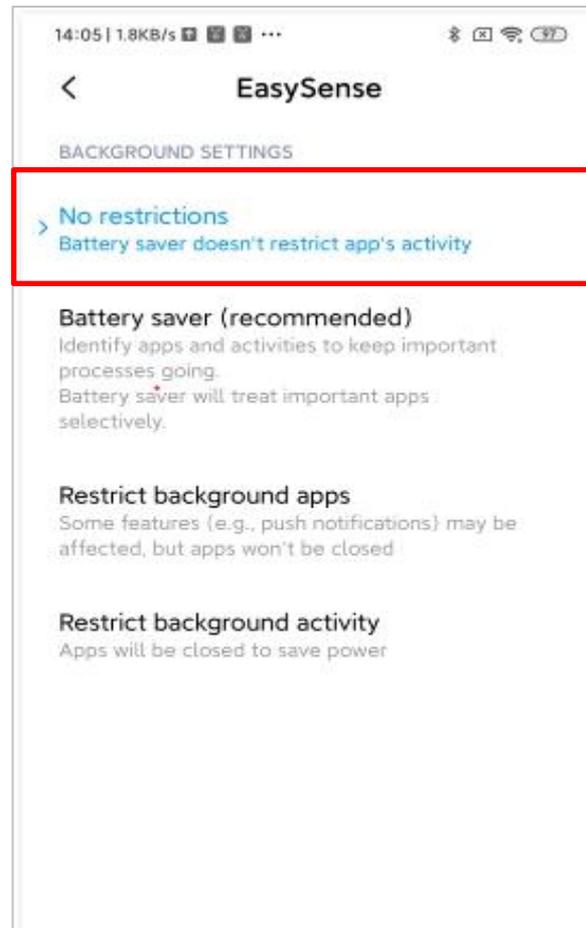

6. Autostart Setting

- Turn on **Autostart** for Medtrum App so that the App is always open and running in the background.

Settings→**Apps**→**Manage apps**→**MedtrumApp**→**Autostart**

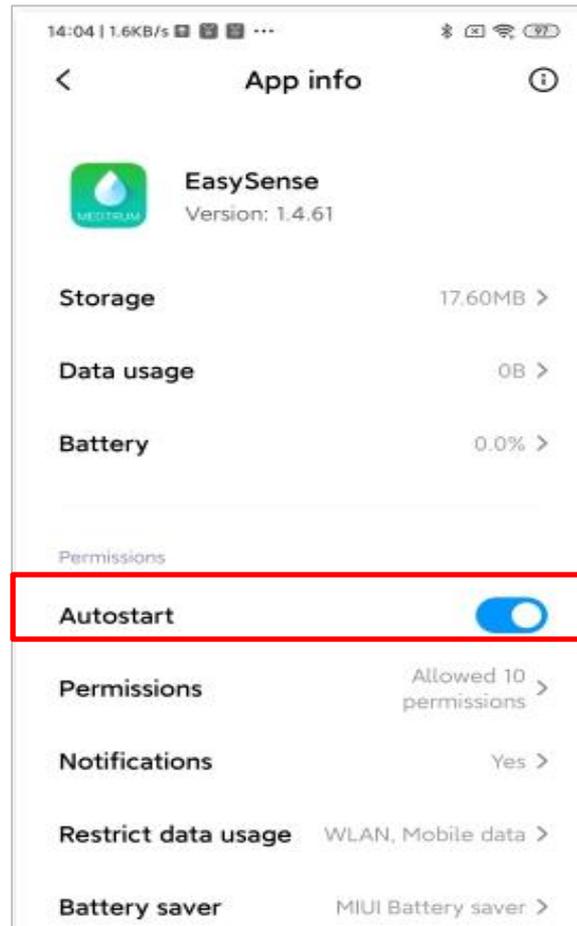

7. Lock Setting

- Tap the lock icon to lock Medtrum App in the background.

